



Hemoglobin Microplate Assay Kit

User Manual

Catalog # ASK1113

Detection and Quantification of Hemoglobin (Hb) Content in Serum Samples.

For research use only. Not for diagnostic or therapeutic procedures.

Bioworld Technology, Inc. (USA)

Email: info@bioworld.com

Web: www.bioworld.com

Bioworld technology, co. Ltd. (China)

Email: info@biogot.com

Web: www.biogot.com



I. INTRODUCTION.....2

II. KIT COMPONENTS.....3

III. MATERIALS REQUIRED BUT NOT PROVIDED.....3

IV. SAMPLE PREPARATION.....4

V. ASSAY PROCEDURE.....5

VI. CALCULATION.....6

VII. TYPICAL DATA.....7



I. INTRODUCTION

Hemoglobin, also spelled haemoglobin and abbreviated Hb, is the iron-containing oxygen-transport metalloprotein in the red blood cells of the blood in vertebrates and other animals. In mammals the protein makes up about 97% of the red cell's dry content, and around 35% of the total content (including water). Hemoglobin transports oxygen from the lungs or gills to the rest of the body, such as to the muscles, where it releases its load of oxygen. Hemoglobin also has a variety of other gas-transport and effect-modulation duties, which vary from species to species, and may be quite diverse in invertebrates.

The Hemoglobin can react with O-tolidine. The products can be measured at a colorimetric readout at 435 nm.



II. KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Dye Reagent	Powder x 1	4 °C
Dye Reagent Diluent	5 ml x 1	4 °C
Reaction Buffer	10 ml x 1	4 °C
Standard	Powder x 1	4 °C
Technical Manual	1 Manual	

Note:

Dye Reagent: add 5 ml Dye Reagent Diluent to dissolve before use.

Standard: add 1 ml distilled water to dissolve, then add 0.1 ml standard solution into 0.9 ml distilled water, the concentration will be 500 µg/ml.

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 435 nm
2. Distilled water
3. Pipettor
4. Pipette tips
5. Centrifuge
6. Timer



IV. SAMPLE PREPARATION

1. For serum sample

Detect directly.



V. ASSAY PROCEDURE

Add following reagents into the microplate:

Reagent	Blank	Standard	Sample
Reaction Buffer	100 µl	100 µl	100 µl
Distilled water	2 µl	--	--
Standard	--	2 µl	--
Sample	--	--	2 µl
Dye Reagent	50 µl	50 µl	50 µl
Mix, measured at 435 nm and record the absorbance.			



VI. CALCULATION


1. According to the serum sample

$$\begin{aligned} \text{Hb } (\mu\text{g/ml}) &= C_{\text{Standard}} \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \\ &= 500 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \end{aligned}$$

C_{Standard} : the concentration of Standard, 500 $\mu\text{g/ml}$.

VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

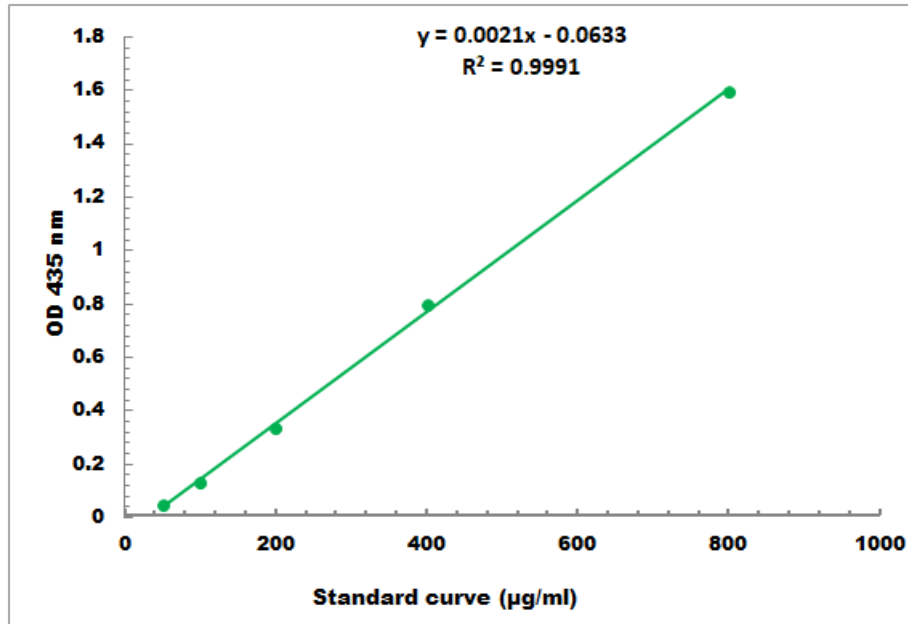

Detection Range: 50 µg/ml - 500 µg/ml